

## Proven over many years:



### Individio Lux

- Order no. 2402 50 tray blanks, blue opaque, upper
- Order no. 2403 50 tray blanks, blue opaque, lower
- Order no. 2408 50 tray blanks, transparent blue, upper
- Order no. 2409 50 tray blanks, transparent blue, lower



### Profibase

- Order no. 2111 25 each of upper and lower baseplates
- Order no. 2112 50 baseplates, upper
- Order no. 2113 50 baseplates, lower



### Individio Light Box

- Order no. 9030 Unit with 4 special bulbs
- Order no. 9031 Replacement bulb

Available from

VC 8-000 2402 GB 094

Contact your dental depot or dental representative and see for yourself the high quality of our products.



## A solid base in the laboratory

A powerful trio:

- Individio Lux
- Profibase
- Individio Light Box

Proven tray blanks and base plates for partial, full and implant prosthetics



P.O. Box 767 · 27457 Cuxhaven · Germany · Tel. +49 (4721) 719-0 · Fax +49 (4721) 719-140 · www.voco.com

creative in research



creative in research

## Light-curing tray blanks and baseplates:

The light-curing tray blanks and baseplates consist of a highly cross-linked resin, which allows a low, uniform plate thickness of 2 mm. Individual colouring and two opacities allow a number of application options. Strongly reduced surface stickiness makes the

trays and base-plates pleasant to handle and to adapt to the model.

Due to its optimum filler content chip removal is the same as with denture acrylics when contouring or trimming with a rotary cutter.

## ● Individuo Lux

### Indications


- Fabrication of custom impression trays for taking functional impressions in partial, full and implant prosthetics
- Bite blocks
- Centric registrations e.g. in Gerber prosthetics
- Customized anterior guide plates in articulators

### Advantages

- Dimensionally stable and fracture-resistant
- Accurate, non-sticky, quick adaptation
- Markings of the model analysis remain visible with transparent material
- No mixing to a dough, ensuring consistent material quality
- Minimum shrinkage
- Light-curing with halogen and UV-A light
- Uniform plate thickness (approx. 2 mm), no reduction required
- Peppermint taste increases patient acceptance
- Low polymerization heat does not melt block out wax

• Even after polymerization Individuo Lux and Profibase have non-sticky, smooth surfaces.

• VOCO tray blanks and base-plates can be easily trimmed with cutters and abrasive rotary instruments.

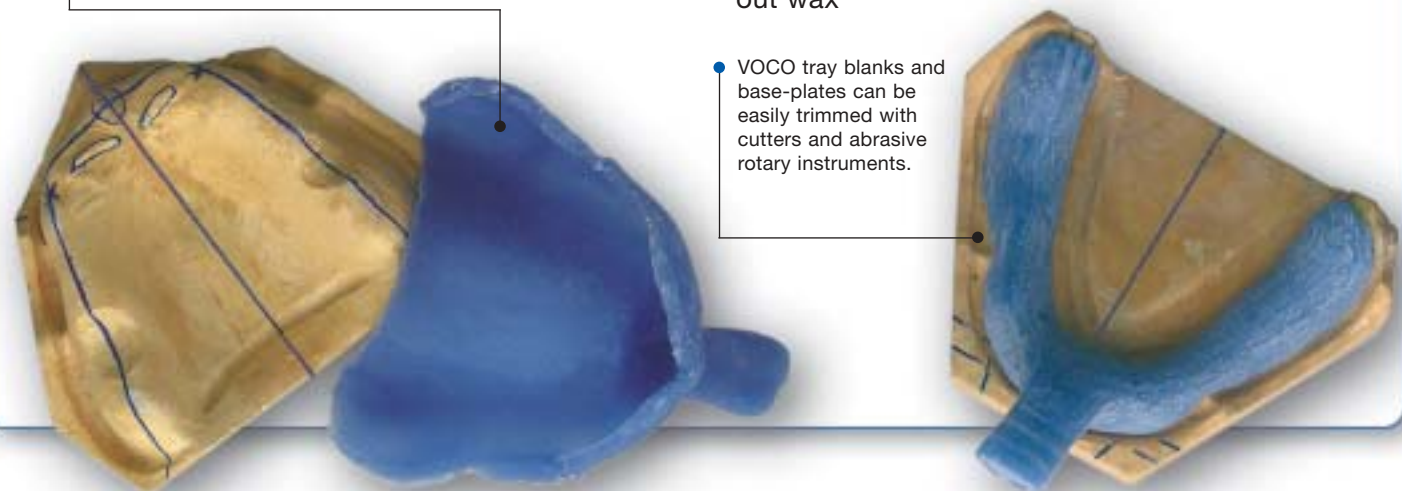

## ● Profibase

### Indications

- Fabrication of custom base plates for setting up partial and full dentures
- Gingival-coloured bite blocks
- Fabrication of surgical stents in implant prosthetics

### Advantages

- Saves time when setting up in wax, no time-consuming adaptation of the wax base
- No deformation of the wax set-up during try-in
- Uniform plate thickness (approx. 2 mm) is the same as the finished denture, avoids time-consuming trimming of the denture base
- Initial assessment of the red and white aesthetics with partial dentures and implant restorations



• Ideal basis for a perfect prosthetic restoration: Profibase

## ● Individuo Light Box

### Indications

- The Individuo Light Box is a robust, easy-to-use halogen light-curing unit for polymerizing Individuo Lux and Profibase tray materials and other light-curing dental materials.

### Advantages

- Easily operated using 2 time programmes
- High polymerization capacity provided by specially designed polymerization chamber
- Spacious sliding tray for up to 4 models is efficient and saves time

### Technical data:

Voltage: 220 V/ 50 Hz

Polymerization light: homogeneous in the blue range, wavelength range 400 - 500 nm

Output: 36 watts (4 x 9 watt bulbs)

W x H x D: 1500 mm X 340 mm x 280 mm